AIR SCALER

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MICRODONT

INDEX | INDICE | ÍNDICE

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Congratulations! You have just purchased a product developed using the latest high-tech concepts available on the market. This equipment has been designed following strict quality standards to provide maximum performance, quality and user safety.

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TECHNICAL NAME:
Dental hand instrument.

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ATTENTION | ATENCIÓN | ATENÇÃO

Read all the information contained in this instruction manual carefully before using the equipment.

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- Cutting-edge elliptical movement, effective in removing dental calculus.
- Its stylized design, anatomical handle and low weight reduce fatigue during long treatments.
- The cooling system provides a mild and constant temperature.
- Easy to install.

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The product should only be used by qualified professionals.


Installation of the Microdont Scaler is very simple. First check that the connection of your equipment is the same as that of the Microdont Scaler Sonic Cavitator (Borden connection), adjust the fitting and tighten it slightly. Once these operations have been carried out, your cavitation instrument will be installed. Next, choose which insert you would like to use and screw it in clockwise by hand. Once this is done, tighten the insert with the appropriate wrench.

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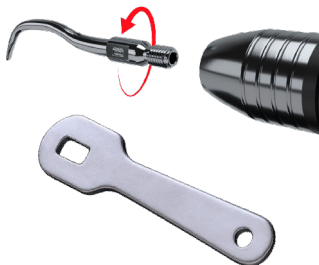ATTENTION:

Excessive effort during the periodontal procedure may cause damage to the tooth.

During use, make sure that the cooling spray is working properly in order to avoid overheating the dental tissues, which could cause damage to the tooth and soft tissues.

Before use, make sure that the inserts are properly secured.

The extractor key is designed to protect the inserts when changing them. To adjust the air pressure and spray jet, proceed as follows: first close the water outlet, press the air drive pedal and adjust the air pressure to 50-60 PSI or 3.4 - 4 BARS. Once the air pressure is set, operate the Microdont Scaler without water until a thin layer of spray forms on the tip.

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WARRANTY | GARANTIA | GARANTIA

The Microdont Scaler Sonic Cavitator is guaranteed for 12 months from the date of purchase against manufacturing defects. The warranty does not apply to normal wear and tear, damage caused by improper use, insufficient maintenance, connection errors, falls or knocks. Keep this certificate very carefully.

In the event of a warranty claim, contact Microdont with the invoice and warranty certificate in hand.

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SPECIAL STORAGE/CONSERVATION/HANDLING CONDITIONS |

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The Microdont Scaler Sonic Cavitator is individually packaged in a case, duly identified and then packed in boxes, thus allowing proper storage and transportation in order to prevent the instrument from being damaged during transportation.

The packaging containing the product must be transported with care, stored and kept in a ventilated, dry place, protected from light and away from the elements.

The products should be kept in their original packaging until they are used.

Handling should only be carried out by accredited persons.

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In order to avoid microbial contamination, this equipment must be sterilized after each use.

- This instrument is fully autoclavable.
- Before autoclaving the Microdont Scaler, all instruments soiled with blood, saliva or other substances must be washed/cleaned.
- Minimum service life: 250 cycles.

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IMPORTANT:

Clean the Microdont Scaler Sonic Cavitator before autoclaving.

The Microdont Scaler Sonic Cavitator must be disposed of in accordance with the legislation in force in the country regarding the disposal of medical/dental materials. Where there is no applicable legislation, it should be disposed of as medical device waste.

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ATTENTION:

In order to guarantee the perfect functioning of the instrument, the air must arrive dry and free of impurities, otherwise it will be considered improper use and any defects caused will not be covered by the warranty.

For best performance, use the correct pressure.

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In case of difficulty, before sending your instrument for technical assistance, check:

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ENGLISH

PROBLEM	CAUSE	SOLUTION
No vibration	Air is not reaching the instrument/Equipment is stuck.	Check the air pressure from the compressor and the unit. Ensure it is within the recommended range.
Low vibration	Low air pressure/Loose tip/Deformed tip	Check the air pressure from the unit and ensure it is within the recommended range/ Tighten the tip using the wrench/Replace the tip with a new one
No water flow	Water passage blocked	For insert blockages: Direct a stream of air and water into the insert orifice/For instrument connection blockages: Direct a stream of air and water into the instrument connection
Water leaking from the connector	Defective gasket	Replace the rubber ring on the coupling

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