ENGLISH

Indication:

The diamond tips for periodontics -Perio Pro- are rotating instruments with AC fitting, which, when connected to a handpiece with low rotation (micromotor + contra-angle), are specifically intended for some periodontal procedures, such as odontoplasty and root planing/root plasty. The kit consists of different formats that are available in granulations: medium, fine and extra-fine, totaling 12 pieces.

Intended user:

Skilled dental professionals — especially periodontists and general practitioners — who perform periodontal procedures involving tooth surface remodeling, root planing, and odontoplasty.

Target population:

Patients undergoing periodontal interventions that require remodeling of tooth surfaces, root planing, or odontoplasty procedures. It includes individuals with periodontal disease, anatomical changes that require regularization, or cases in which root refinement contributes to periodontal control and better therapeutic response.

Rationale for action:

The diamond tips of the *Perio Pro* kit act by controlled abrasion of the root surface. Its diamond coating, associated with rotation at low speed, allows the precise removal of irregularities, adhered deposits and contaminated surface areas. This process promotes the smoothing and regularization of the treated structures, favoring the reduction of biofilm retention.

Contraindication:

- Do not use in patients with severe tooth mobility.
- Do not use at high speed.
- Cleaning agents with chlorine or chloride as the active ingredient are corrosive and should not be used.
- Do not use cold sterilization methods (chemical process). These chemical agents can corrode the product.


Instructions for use:

- Insert the tip into a low-speed handpiece (micromotor + contra-angle), which meets the applicable requirements and standards for the product.
- Before activating the handpiece, check that the tip is properly fitted and attached to the instrument.
- The instrument must not exceed 6,000 RPM, pressure must not be applied during use, and coolant must be engaged at all times.
- At the end of the procedure, proceed with the cleaning and sterilization of the tip.

Form of presentation

- Non-sterile product;
- Kit with bur block + 12 diamond tips;

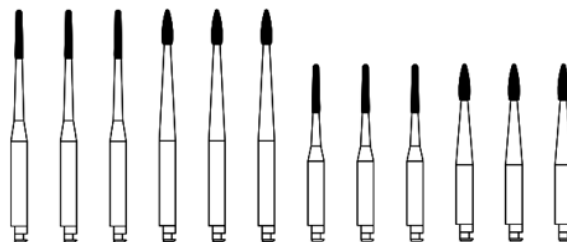Cleaning procedures (Kit- Perio Pro):

- Clean and sterilize the diamond tips before each use according to the cleaning and sterilization item.
- It is recommended that cleaning be done immediately after the end of the procedure, to avoid drying out contaminants, which can make cleaning difficult.
- After cleaning, check the tips to ensure that any traces of contamination have been removed. Repeat cleaning if necessary.
- Cleaning of diamond surgical tips can be done in two ways: Manual or Ultrasonic.

Manual Cleaning Procedure:

If manual cleaning is the only option available, proceed as follows.

- Use an area specifically reserved for cleaning;
- Rinse the diamond tips under running water and keep them immersed in a bath of neutral pH cleaning solution (enzymatic detergent) following the agent manufacturer's instructions.
- After soaking and keeping it immersed, brush well away from the body using the neutral cleaning agent (enzymatic detergent). Care should be taken to avoid spreading contaminants by spraying during the brushing process. If necessary, it is possible to use specific brushes for cleaning dental drills/tips.
- Rinse the diamond tips thoroughly under running water until they are visibly clean.
- Dry the device with clean compressed air or paper towels.

Ultrasonic Cleaning Procedure:

For ultrasonic cleaning, it is recommended to proceed as follows:

- Prepare a pH-neutral cleaning solution (enzymatic detergent) following the agent manufacturer's instructions for correct concentration, exposure time, temperature, and water quality.
- Place the diamond tips for periodontics in the appropriate location for the ultrasonic equipment, following the manufacturer's instructions for use.
- Completely immerse the device in the cleaning solution and program the equipment for at least fifteen (15) minutes.
- Rinse the diamond surgical tips thoroughly under running water until they are visibly clean.
- Dry the device with clean compressed air or paper towels.

Sterilization requirements for diamond tips:

- It is recommended to sterilize the Perio Pro Kit in an autoclave with saturated steam under pressure at 121°C for 20 minutes, or to follow the sterilization parameters according to the relevant policies/regulations of the country where the device is being used.
- It is the professional's sole responsibility to keep the autoclave equipment up to date to ensure the correct sterilization of the product.

USE LIMITATIONS: The recommended shelf life of the product is 5 uses (5 reuses), however it is of paramount importance that before each use a check is carried out, that is, if the dental diamond tips are in perfect condition for use, otherwise they can overheat the dental tissues and lead to an unwanted pulp response triggering a painful sensation in the patient, in addition to not complying with the necessary wear requirement.

Warnings and Precautions:

- Clean and sterilize dental diamond tips before each use according to cleaning and sterilization items
- For user and patient safety, the dental diamond tip should only be used by qualified dental professionals.
- Before use, check that the dental diamond tips are in perfect condition for use, especially in relation to the diamond area and integrity of the rod
- Before use, make sure that the rod of the dental diamond tips is fully seated in the low-rotation handpiece;
- Never force dental diamond tips into the low-speed handpiece. In case of difficulty in fitting, check the turbine and the dental diamond tips.
- Respect the working speed (RPM Max) defined on the product label, otherwise you may generate undesirable heat.
- Chlorine or chloride cleaning agents as an active ingredient are corrosive to stainless steel and should not be used.
- Prolonged immersion of dental diamond tips in disinfectant solutions can result in corrosion and should therefore be avoided.
- Do not use cold sterilization methods (chemical process) for sterilization of dental diamond tips. These chemical agents can cause corrosion of the products.

Using:

- Always act with intermittent movements;
- Move the dental diamond tips continuously to avoid localized heating and consequent breakage of the tips, causing risks to the patient;
- Use proper cooling methods to ensure proper cooling of the parts used for treatment;
- Use aspiration techniques that do not compromise the cooling of the parts involved in the treatment;
- Eye protection should be worn to protect against ejected particles.
- Used dental diamond tips should be considered contaminated instruments and appropriate handling precautions are required. Wear personal protective equipment (PPE), including gloves, masks, and safety glasses when handling, as with all other dental instruments

Conservation and storage:

- The product must be stored in an airy, dry place protected from light, in order to maintain ideal storage and transport conditions, as well as its physical and chemical integrity.
- Keep the product in its original packaging until the moment of use. Inspect the packaging before use.
- The carrier must be informed about the contents. The precautions to be taken in transport are demonstrated on the packaging, thus ensuring the protection of the product from shipping to delivery to the customer. Packets should not be left to the action of time.
- Transport in a dry, clean place at room temperature.

Discard:

The disposal of the material must follow the relevant local legislation.

Shelf Life:

Not applicable. The product is manufactured from non-perishable materials and supplied non-sterile. When stored under the conditions indicated, it retains its properties indefinitely.

Additional information:

- If necessary, 3R provides instructions for use in printed form, at no additional charge, upon request by email (qualidade@3r.ind.br) or by phone: (+55 11-5525-0590 / +55 0800-773-5525), and we will get back to you within a maximum of 7 business days from the request.
- Always check the website: www.microdont.com.br, the most up-to-date version of the instructions for use corresponding to the product.

Composition:

AISI 303 stainless steel
Diamond

Symbology;

Symbol	Symbol Title	Description
	Validity	Indicates the date after which the health product should not be used.
	Non-sterile	Indicates a health product that has not undergone a sterilization process.
	Keep dry	Indicates a health product that needs to be protected against moisture.
	Medical device	Indicates that it is a medical device.
	Recyclable packaging	The packaging is recyclable.
	See instructions for use	Indicates the user's need to refer to the instructions for use.
	Date of manufacture	Indicates the date the health product was manufactured.
	Keep out of sunlight	Indicates a health product that needs protection from light sources.
	Manufacturer	Indicates the manufacturer of the health product.
	Lot Code	Indicates the manufacturer's batch code so that the batch or lot can be identified.
	Catalog Number	Indicates the manufacturer's reference number so that the health product can be identified.
	For professional use only	It indicates that the medical device must be used exclusively by a qualified health professional.
Rx Only	Prescription only	Sale and use restricted to licensed professionals or by prescription, in accordance with applicable legislation.
	Unique Device Identification (UDI)	Indicates that the adjacent code corresponds to the Unique Device Identifier (UDI).
	QR Code – Access to the electronic IFU	QR Code to access the electronic version of the product's Instructions for Use (eIFU).

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