ENGLISH**Indication:**

Polishing discs are indicated for finishing and polishing restorations in composite resin and glass ionomer, and are designed to provide a smooth surface. For this, they should be used starting with the coarsest grain and moving on to the finest. They act as auxiliaries in the finishing and polishing stages, which must be complemented with the use of other polishers.

Intended user:

Duly qualified dental professionals.

Target population:

The product is intended for use in patients who have undergone restorative procedures, with the aim of providing the proper finishing and polishing of the restoration.

Fundamentals of its operation and action:

The polishing discs act by means of controlled abrasion on the surface of the restorative material. The abrasive particles present on the surface of the disc promote the gradual removal of irregularities and grooves formed during the adjustment or finishing of the restoration. The sequential use of different grains, from the coarsest to the finest, allows the progressive reduction of surface roughness, contributing to obtaining a smoother surface with better final gloss.

Applications:

Finishing and polishing of restorations made in composite resin or glass ionomer.

Contraindications:

- The product should not be used in patients with a history of hypersensitivity to any components present in its composition.
- Do not use in restorations made with materials other than composite resin or glass ionomer.
- Do not use the product for purposes other than those for which it was developed.


Instructions for Use:

- The Polidont Polishing Discs come with a metal rod to be fitted and used in the low rotation motor.
- Fit the discs onto the metal rod.
- Fit the mandrel shank to the low-rotation handpiece (counter angle).
- Make sure the metal rod and disc are properly seated in the low-rpm motor.
- Check that the refrigeration is working properly.
- Start the motor and wear/polish the desired region.

The sequence of the procedure for the use of the discs should be in the decreasing direction from higher grain (extra coarse) to the lowest grain size (ultrafine).

POLIDONT
POLIDONT
POLIDONT

REV.

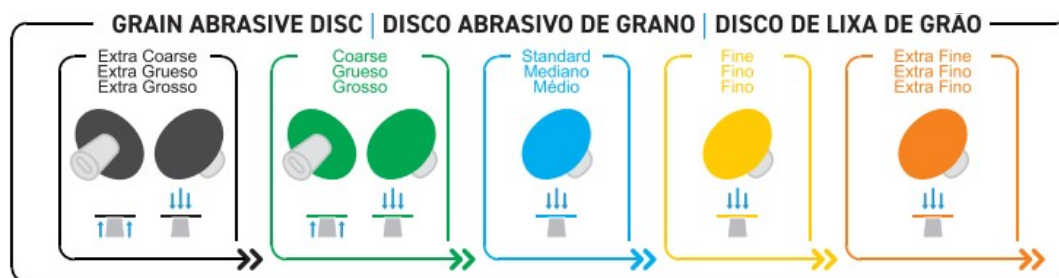BLACK – EXTRA COARSE GRAIN:

Indicated for the initial finish of the restoration, with the purpose of defining the shape and contour of the tooth and removing excess restorative material.

Due to the greater rigidity of the disc, it is recommended to use it for the definition of buccal and incisal embrasures. Because it has extra coarse grain, it promotes greater abrasion of the restorative surface, which can generate deeper grooves. It is recommended to carry out irrigation or washing of the surface after use, to remove aluminum oxide particles and prevent excessive wear.

The discs can have an active inner or outer face:

- Active outer face: suitable for finishing smooth faces.
- Active inner face: indicated for finishing incisal and/or interproximal faces.

GREEN – COARSE GRAIN: Also indicated to be used in the finishing phase of the restoration. Its surface with coarse abrasive, when applied to the restored surface, promotes the reduction and flattening of the irregularities left by the previous disc.

It can have an active internal or external face, with the same indications for use described above:

- Active outer face: smooth face finish.
- Active inner face: finishing of incisal and/or interproximal faces.

BLUE – MEDIUM GRAIN: Indicated for the final finish of the restoration, being used to reduce or eliminate remaining surface grooves from the previous finishing steps. During this stage, it is recommended to use an intermittent water jet in order to favor the removal of abrasive aluminum oxide particles and avoid heating of the restored surface.

YELLOW – FINE: Suitable for the initial polishing of the restoration. It should be used with a slightly moistened tooth surface. Lubricant can be used during the procedure. The use of continuous waterjet is not necessary due to the small diameter of the abrasive particles.

ORANGE – EXTRA FINE: Indicated for the final polishing of the restoration, promoting greater surface smoothness and final shine of the restorative material.

Form of presentation:

- Non-sterile product.
- Single use.
- Black: Extra Coarse / Green: Coarse / Blue: Medium / Yellow: Fine / Orange: Extra Fine.
- Polidont Introductory Kit: 28 discs + 1 Mandrel
- Polidont Introductory Mini Kit: 7 discs + Mandrel
- Polidont Introductory Mini Kit Refill: 7 discs

Cleaning Procedure:

- The discs are single-use; therefore, they should not be cleaned or reused.
- Only the mandrel can be cleaned and reused.

Sterilization Procedure:

- The discs are supplied non-sterile and should not be sterilized. They should be disposed of after use, according to local regulations for dental waste.
- Only the mandrel can be sterilized and reused

Warnings and Precautions:

- To avoid fractures and distortion of the instrument, do not apply excessive force when inserting the disc into the low-speed handpiece.
- Before use, inspect if the disc is in perfect condition for use, especially in relation to the area with granulation and the integrity of the rod;
- Before use, make sure the disc rod is perfectly fitted to the contra-angle;
- Never force the disc into the counter-angle (low rotation). In case of difficulty in fitting, check the contra-angle (low rotation) and the rod of the Disc.
- It is recommended that the use of the product is always accompanied by absolute isolation of the operative field.
- Do not use the disc as levers and/or at incorrect angles;
- Contact with soft tissues (gums, mucosa) can cause local lesions.
- Cleaning agents with chlorine or chloride as the active ingredient are corrosive to stainless steel and should not be used.
- Prolonged immersion of the Disc in disinfectant solutions may result in corrosion and should therefore be avoided.
- Protective eyewear should be worn to protect against particles ejected on the patient and the professional.
- Wear PPE, including gloves, masks, safety goggles, and long-sleeved lab coats when handling the polishing disc, as is done with all other dental instruments. Wear protective glasses on the patient as well.
- In case of returning the product, make sure that it has been effectively cleaned and properly packaged for return.
- BATCH NUMBER: It is important to keep the batch number, identified in the product title, to be informed of any occurrences with the manufacturer.

Conservation and storage:

- Keep the product in its original packaging until the moment of use. Inspect the packaging before use.
- The product must be stored in an airy, dry place protected from light, in order to maintain ideal storage and transport conditions, as well as its physical and chemical integrity.
- The carrier must be informed of the contents. The precautions to be taken in transport are demonstrated on the packaging, thus ensuring the protection of the instrument from shipping to delivery to the customer. Packets should not be left to the action of time.

Discard:

POLIDONT
POLIDONT
POLIDONT

REV.

- The disposal of the material must follow the relevant local legislation.

Expiration date:

- Not applicable. The product is manufactured from non-perishable materials and supplied non-sterile. When stored under the conditions indicated, it retains its properties indefinitely.

Shelf life:

- Not applicable. The product is for single use only and should be discarded after use.

Additional information:

- If necessary, 3R provides instructions for use in printed form, at no additional charge, upon request by email (administrativo@3r.ind.br) or by phone: (+55 11-5525-0590 / +55 0800-773-5525), and we will get back to you no later than 7 business days from the request.
- Always check the website: www.microdont.com.br, the most up-to-date version of the instructions for use corresponding to the product.

Composition:

Active Part - Polyester and Aluminum Oxide

Mandrel – AISI 304 stainless steel

Symbology:

Symbol	Symbol Title	Description
	Validity	Indicates the date after which the health product should not be used.
	Non-sterile	Indicates a health product that has not undergone a sterilization process.
	Do not reuse	Indicates an intended health product for use or for use on a single patient during a single procedure.
	Keep dry	Indicates a health product that needs to be protected against moisture.
	Medical device	Indicates that it is a medical device.
	Recyclable packaging	The packaging is recyclable.
	See instructions for use	Indicates the user's need to refer to the instructions for use.
	Date of manufacture	Indicates the date the health product was manufactured.
	Keep out of sunlight	Indicates a health product that needs protection from light sources.
	Manufacturer	Indicates the manufacturer of the health product.
	Lot Code	Indicates the manufacturer's batch code so that the batch or batch can be identified.
	Catalog Number	Indicates the manufacturer's reference number so that the health product can be identified.
	For professional use only	It indicates that the medical device must be used exclusively by a qualified health professional.
Rx Only	Prescription only	Sale and use restricted to licensed professionals or by prescription, in accordance with applicable legislation.
	Unique Device Identification	Indicates that the adjacent code corresponds to the Unique Device Identifier (UDI).

POLIDONT
POLIDONT
POLIDONT

REV.

Symbol	Symbol Title	Description
	(UDI)	
	QR Code – Access to the electronic IFU	QR Code to access the electronic version of the product's Instructions for Use (eIFU).

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